



KERN FXN-N

Robust, space-saving bench scale, protected by stainless steel and IP68 rating – now with EC type approval [M]

Features

- Suitable for the ever-increasing hygienic requirements in the food industry
- Your support in a **HACCP**-compliant quality system
- Ideal for the robust industrial applications
- **protection against dust and water splashes IP68**, ideal for harsh ambient conditions. Submersion in water actually possible
- **Stainless steel design** of the housing and weighing plate. Its smooth surfaces make it simple to clean
- **Very fast display**: steady weight values within 2 s
- **High mobility**: thanks to battery operation, compact, lightweight construction, it is suitable for the use in several locations (kitchen, sales office, cafeteria, Food industry-Laboratory etc.)

• **1 Recessed grips on the underside**

- of the scale for easy transportation
- **Weighing with tolerance range (check-weighing)**: a visual and audible signal helps with portioning, dispensing or grading

Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions weighing surface, stainless steel, WxD 236x195 mm
- Overall dimensions WxDxH 240x280x120 mm
- Ready for use: Batteries included, 4x1.5 V Size D, operating time up to 200 h with backlight, AUTO-OFF function to preserve the battery
- Net weight approx. 3,2 kg
- Permissible ambient temperature 0 °C/40 °C

Accessories

- **2 Tare pan made of stainless steel**, ideal for weighing loose small parts, fruit, vegetables etc., WxDxH 370x240x20 mm, KERN RFS-A02
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD



OPTION


























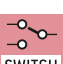




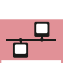










Model	Weighing range [Max] kg	Readout [d] g	Reproducibility g	Linearity g	Option	
					DAkkS Calibr. Certificate	
KERN					DKD KERN	
FXN 3K-4N	3	0,5	0,5	± 1	963-127	
FXN 6K-3N	6	1	1	± 2	963-128	
FXN 10K-3N	15	2	2	± 4	963-128	
FXN 30K-3N	30	5	5	± 10	963-128	
FXN 3K-3*	3	1	1	± 2	963-127	↓
FXN 6K-3*	6	2	2	± 4	963-128	↓
FXN 10K-3*	15	5	5	± 10	963-128	↓
FXN 30K-2*	30	10	10	± 20	963-128	↓

! *ONLY WHILE STOCKS LAST !

↓ Price reduction

KERN Pictograms:

 Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).	 Piece counting: Reference quantities selectable. Display can be switched from piece to weight.	 Rechargeable battery pack: Rechargeable set.
 Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required.	 Recipe level A: Separate memory for the weight of the tare container and the recipe ingredients (net total).	 Universal mains adapter: with universal input and optional input socket adapters for A) EU, GB B) EU, GB, CH, USA C) EU, GB, CH, USA, AUS
 Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	 Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display.	 Mains adapter: 230V/50Hz in standard version for EU. On request GB, USA or AUS version available.
 Alibi memory: Electronic archiving of weighing results, complying with the 2014/31/EU standard.	 Recipe level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, multiplier function, adjustment of recipe when dosages are exceeded or barcode recognition.	 Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request.
 Data interface RS-232: To connect the balance to a printer, PC or network.	 Totalising level A: The weights of similar items can be added together and the total can be printed out.	 Weighing principle: Strain gauge Electrical resistor on an elastic deforming body.
 RS-485 data interface: To connect the balance to a printer, PC or other peripherals. High tolerance against electromagnetic disturbance.	 Percentage determination: Determining the deviation in % from the target value (100 %).	 Weighing principle: Tuning fork A resonating body is electromagnetically excited, causing it to oscillate.
 USB data interface: To connect the balance to a printer, PC or other peripherals.	 Weighing units: Can be switched to e.g. non-metric units at the touch of a key. See balance model. Please refer to KERN's website for more details.	 Weighing principle: Electromagnetic force compensation Coil inside a permanent magnet. For the most accurate weighings.
 Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals.	 Weighing with tolerance range: Upper and lower limiting values can be programmed individually for e.g. dosing, sorting and portioning.	 Weighing principle: Single cell technology Advanced version of the force compensation principle with the highest level of precision.
 WLAN data interface: To transfer data from the balance to a printer, PC or other peripherals.	 Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value.	 Verification possible: The time required for verification is specified in the pictogram.
 Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.	 Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.	 DAkKS calibration possible (DKD): The time required for DAkKS calibration is shown in days in the pictogram.
 Interface for second balance: For direct connection of a second balance.	 ATEX explosion protection: Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.	 Package shipment: The time required for internal shipping preparations is shown in days in the pictogram.
 Network interface: For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.	 Stainless steel: The balance is protected against corrosion.	 Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram.
 Wireless data transfer: between the weighing unit and the evaluation unit using an integrated radio module.	 Suspended weighing: Load support with hook on the underside of the balance.	 Warranty: The warranty period is shown in the pictogram.
 GLP/ISO log: The balance displays the weight, date and time, regardless of a printer connection.	 Battery operation: Ready for battery operation. The battery type is specified for each device.	
 GLP/ISO log: With weight, date and time. Only with KERN printers.		

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAkKS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkKS calibration laboratory today is one of the most modern and best-equipped DAkKS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAkKS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAkKS calibration of balances with a maximum load of up to 50 t
- DAkKS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAkKS calibration certificates in the following languages DE, GB, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Your KERN specialist dealer: